

醫卡 funcare®



倍熱代言人  
朱可自

# Massage Health maintenance

funcare® digital tens



## The function of each part

### LCD display window

The massager's LCD display screen displays current treatment program, intensity and duration usage.

### On/Up key

Press this button to switch on the massager. Press once to increase massage intensity. The highest intensity is 15.

### Conductive pads

Please use the jelly-like "sticky patch" with the conductor patches to ensure that the patches stick to your body without falling off.


### The down/off key

Press this button to decrease massage intensity. Press and hold to switch off the massager.

### Program mode selection

Press this button to choose from 10 different massage treatments.

Battery cover

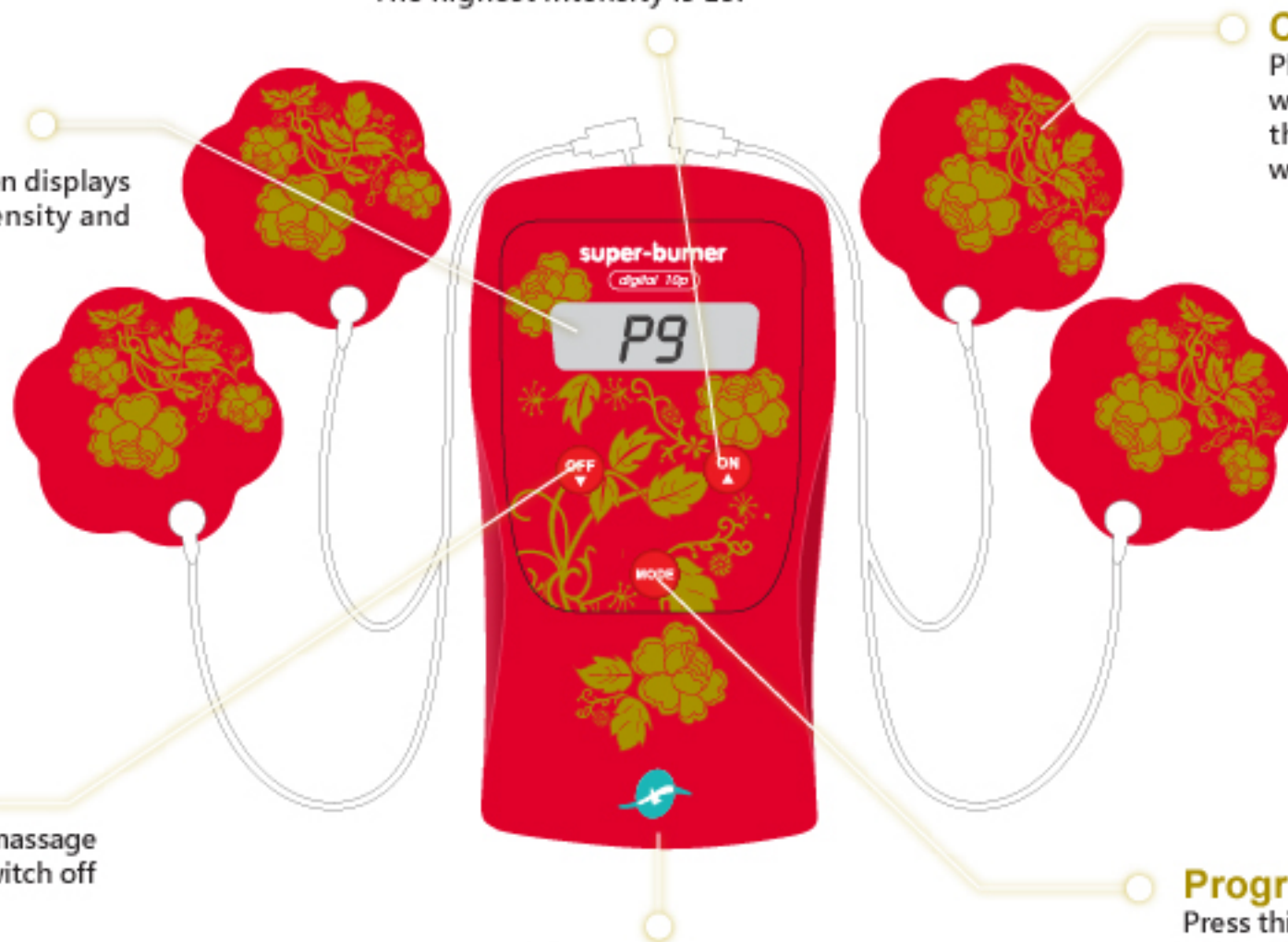

## Installation of batteries


〈 Must use one 9V battery 〉

\* Confirm the main unit is turned off before changing the battery.

Front




Back

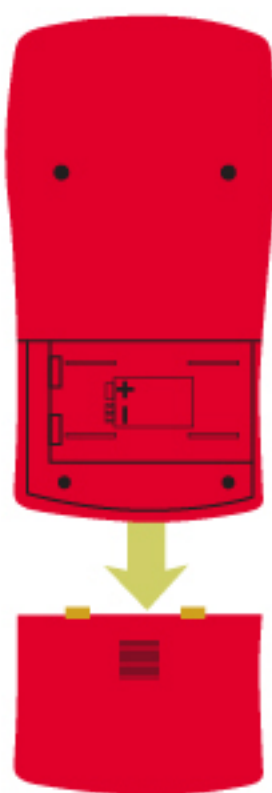

2

Prepare one 9V battery and put the battery into the case according to the instructions on the battery case.

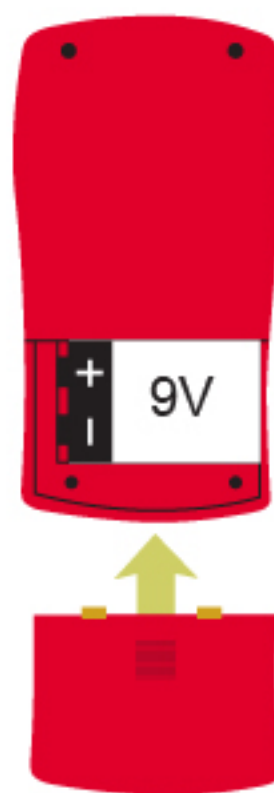

1

Remove the battery case from the back of the massager.

3

Put the battery back to the original position.

## Installation of conductive pads



**1** Insert one side of the conductive cable to the dent of the main unit.



**2**

Connect the other side of the cable to 2 conductive pads to enable the user to easily operate it. The pads can be attached to 4 parts for use concurrently.

