



*Like the sensation from
professional masseurs*

Relaxes and relieves stress

Funcare Digital Tens(mini)



松柏復健科診所 總院長

蘇宗柏 醫師

The function of each part

LCD display window

The massager's LCD display screen displays current treatment program, intensity and duration usage.

Program mode selection

Press this button to choose from 10 different massage treatments.

On/Up key

Press this button to switch on the massager. Press once to increase massage intensity. The highest intensity is 15.



The down/off key

Press this button to decrease massage intensity. Press and hold to switch off the massager.

Battery cover

Conductive pads

Please use the jelly-like "sticky patch" with the conductor patches to ensure that the patches stick to your body without falling off.

Installation of batteries


(Please only use 3V CR2025 lithium batteries)

* Please ensure that the massager is switched off before removing or installing batteries.

1. Remove the bottom part of the battery case from the massager.



2. Place the batteries into the case according to the instructions provided. Hold the batteries on the side; make sure that opposite poles of batteries do not come into contact. This ensures that the batteries do not lose their power.

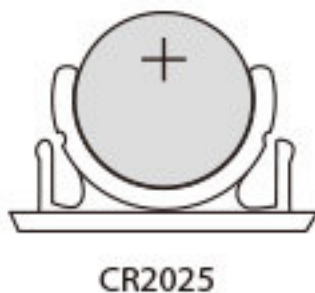


3. Place the battery box back into the massager.



Installation of conductive pads

How to install the conductor patches:

Wireless Method:

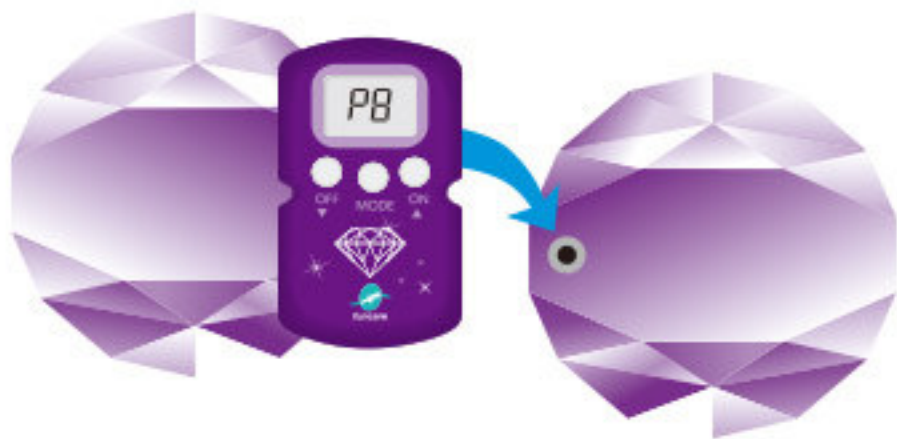

Attach the two conductor patches directly onto back of the massager. The wireless method conceals the patches better and allow for a more convenient and comfortable massage.

Wired Method:

This method requires the use of the*1 to2* wire. Attach the single-end part of the wire to the back of the massager. Attach the double-end part of the wire to the two conductor patches. This method allows user to massage two body parts of a more intense massage experience.

