

Relaxes and relieves stress

Alleviates soreness / Detects aches and pains

Funcare Digital Tens – Aching Muscles Massager (Purple)



蘇宗柏 醫師
松柏復健科診所 總院長



資深醫療法律律師

歐清彬

酸痛專家


廖康吟

The function of each part

LCD display window
The massager's LCD display screen displays current treatment program, intensity and duration usage.

Program mode selection
Press this button to choose from 10 different massage treatments.

Conductive pads
Please use the jelly-like "sticky patch" with the conductor patches to ensure that the patches stick to your body without falling off.

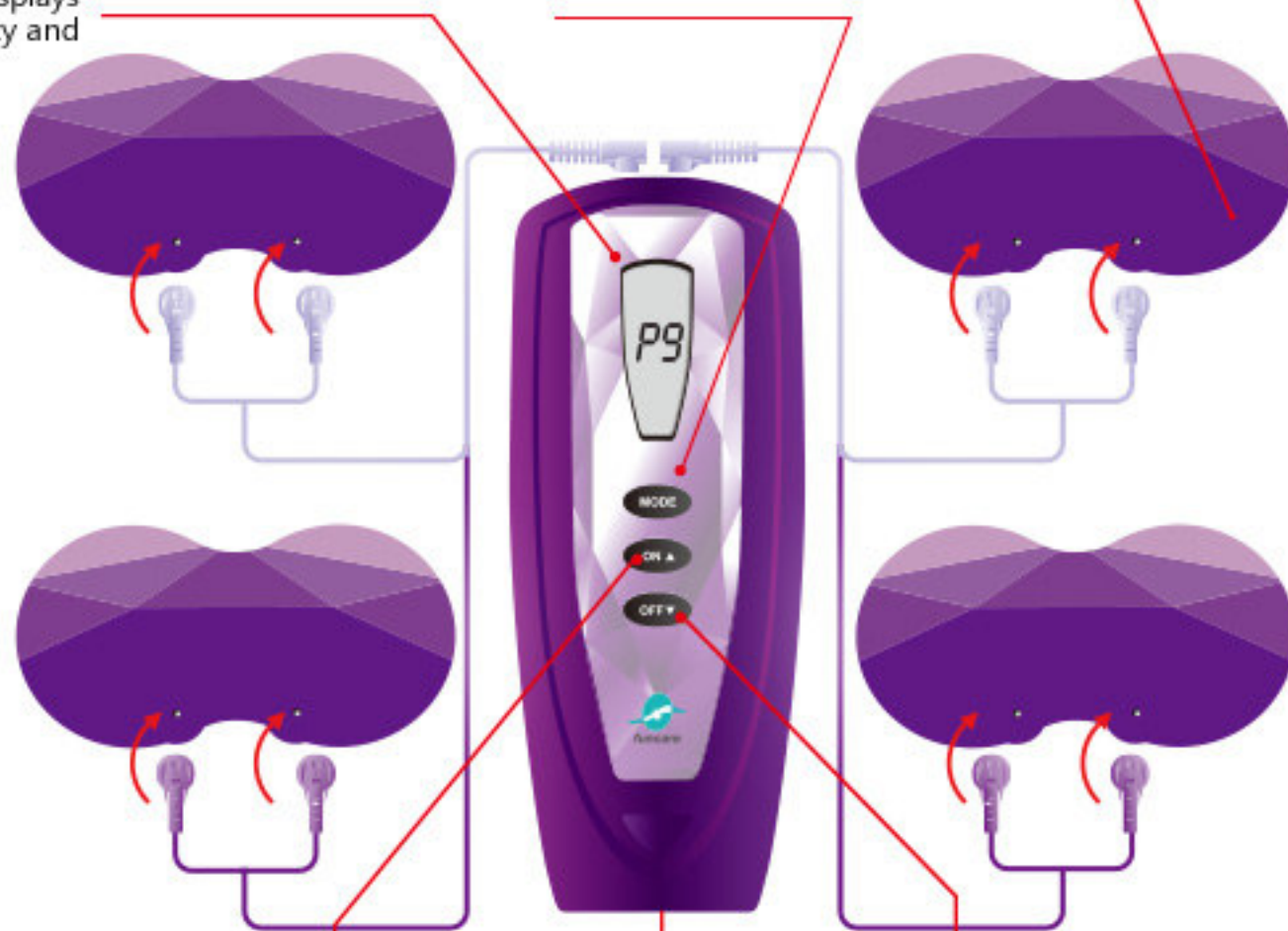

On/Up key
Press this button to switch on the massager.
Press once to increase massage intensity.
The highest intensity is 15.

Battery cover

The down/off key
Press this button to decrease massage intensity.
Press and hold to switch off the massager.

Installation of batteries


〈 Must use one 9V battery 〉

* Confirm the main unit is turned off before changing the battery.

1. Remove the battery case from the back of the massager.



3. Put the battery back to the original position.

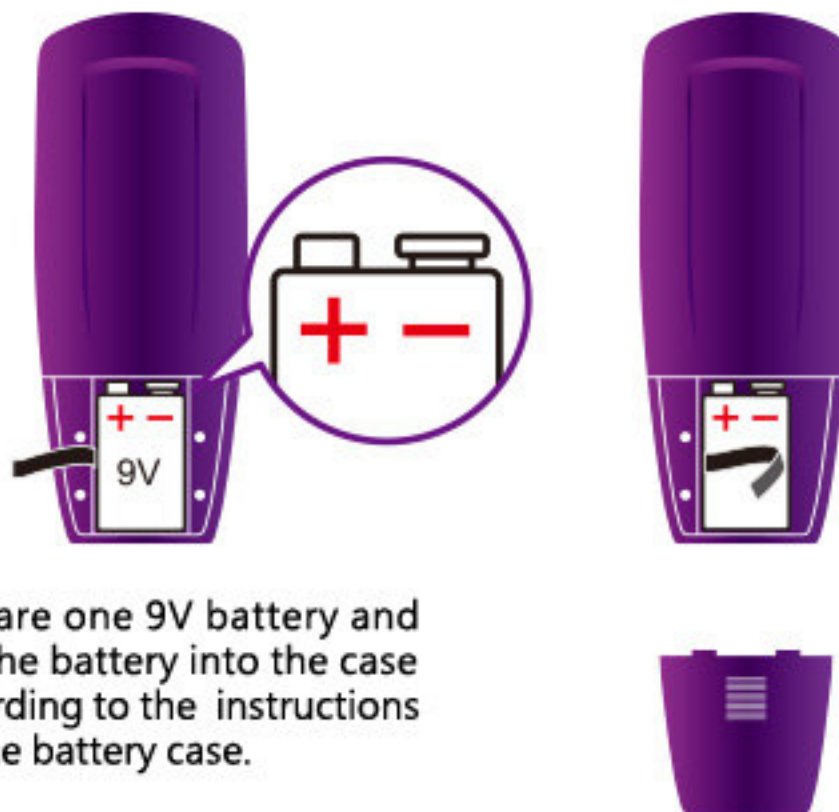

2. Prepare one 9V battery and put the battery into the case according to the instructions on the battery case.

Installation of conductive pads

1. Insert one side of the conductive cable to the dent of the main unit.
2. Connect the other side of the cable to 2 conductive pads to enable the user to easily operate it. The pads can be attached to 4 parts for use concurrently.

